
Baha® Intenso™

User Manual

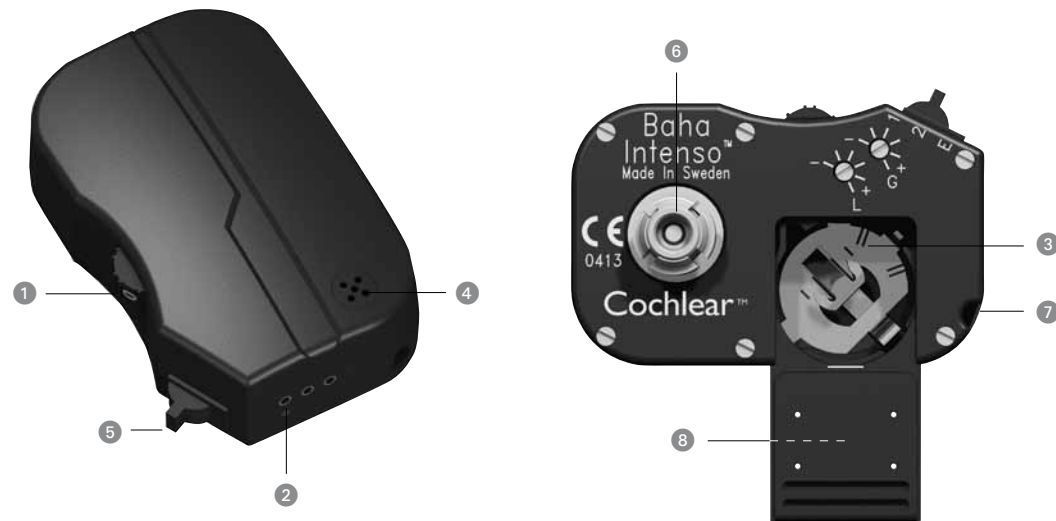Contents

English	4
Français	21
Español	38
Português	55

Baha Intenso™

- | | |
|-----------------------|--|
| ① Volume control | ⑤ Program selector |
| ② Direct audio input | ⑥ Plastic snap connector |
| ③ Battery compartment | ⑦ Attachment point for the safety line |
| ④ Microphone | ⑧ Serial number |

The Baha Intenso™ is a bone conduction sound processor and the most powerful of our head worn devices. The Baha® bone conduction system works by combining a sound processor with an abutment and a small titanium implant placed in the skull behind the ear. The system is based on a process of "osseointegration" through which living bone tissue integrates with titanium. Thus, the titanium implant becomes one with the bone allowing sound to be conducted via the skull bone directly to the cochlea and improving your hearing.

Please read this manual carefully to learn how to use and maintain your sound processor. Be sure to discuss any questions or concerns that you may have regarding hearing or use of this system with your hearing care professional.

We strive to provide the best possible products and services to our clients. Your views and experiences of our products and services are important to us. If you have any comments you would like to share, please contact us.

Address:

Customer Service - Cochlear Americas

13059 E. Peakview Avenue, Centennial, CO 80111, USA
 Toll free (North America) 1800 523 5798
 Tel: +1 303 790 9010 Fax: +1 303 792 9025
 Email: customer@cochlear.com

Customer Service - Cochlear Europe

6 Dashwood Lang Road, Bourne Business Park, Addlestone, Surrey KT15 2HJ, UK
 Tel: +44 1932 87 1500 Fax: +44 1932 87 1526
 Email: info@cochlear.co.uk

Customer Service - Cochlear Asia Pacific

1 University Avenue, Macquarie University, NSW 2109, Australia
 Toll free (Australia) 1800 620 929
 Toll free (New Zealand) 0800 444 819
 Tel: +61 2 9428 6555 Fax: +61 2 9428 6352 or Toll free 1800 005 215
 Email: customerservice@cochlear.com.au


Getting started

Connecting and disconnecting the sound processor

A plastic snap connector is mounted to the sound processor. It is designed to snap onto the abutment and hold the sound processor securely in place.

Connecting ► Hold the sound processor so that the plastic snap connector is accessible and move any hair away from the abutment. Push the sound processor carefully in with a tilting motion and it will attach to the abutment.



Disconnecting ► Tilting the sound processor slightly in any direction will release it from the abutment.

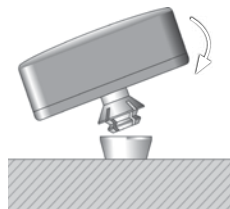Warning! Do not push the sound processor straight in and do not pull it straight out. This puts pressure on the abutment and may cause discomfort.

On/Off

Off ► When the volume control is in the 0 position it is off. You will feel a “click” as it goes into off mode.

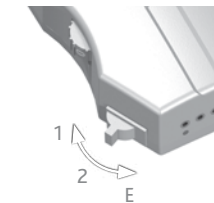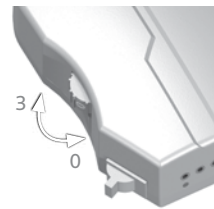On ► For the sound processor to be on, the volume control must be set to a position higher than 0. You will feel a click as the sound processor moves from off mode to on mode.

Volume control ► If you turn the volume control upwards the volume increases, and if you turn it downwards it will decrease.

Program selector

As you change from one program to another sounds signal which program you enter. One beep for program one, two for program two and three beeps for program E.

Program 1 ► Program one is the “normal setting”, suitable for most everyday sound environments.

Program 2 ► For comfort and convenience program two has a higher bass cut setting that some people will find useful for reducing background noise in loud environments.

Program E ► This program is used when connecting electrical accessories such as the audio adapter or telecoil. If no accessory is connected, there will be no sound from the device even if it is on.

Gain and Tone control

The gain and tone controls are to be adjusted **by audiologists only**.

When fitting the Intenso, the audiologist will set the gain and tone controls to suit your needs by increasing or decreasing the volume and/or low frequency outputs.

Low battery warning

The Intenso has a low battery warning that will signal the need to change battery. As the battery is losing power you will hear two longer “beeps”.

Telephones and mobiles

When using telephones, including mobiles, with your Baha Intenso™, you hold the receiver in the normal manner, but close to the microphone inlet of the sound processor instead of close to your ear. It is important that neither you, nor the phone, touches the sound processor as this may cause feedback (“whistling”).

Care and maintenance

Cleaning the abutment area

When the last dressing has been removed and the abutment is exposed hygiene continues to be very important.

Daily care

On a daily basis the area around the abutment should be cleaned to avoid debris build up. This build up around the base of the abutment is not a scab and must be removed. Cleaning this area is most easily done when you bath or take a shower as plenty of warm water and soap on the area will help to soften any crust that may have developed around the base of the abutment. A soft brush is provided together with your sound processor for this purpose.

During the first weeks, after the dressing is removed, cleaning with the soft brush may feel a little rough. If so, use a strip of baby wipe until your skin becomes less sensitive and you can use the brush comfortably. The soft brush can be used with any mild soap or shampoo and plenty of warm water.

When you brush, aim the bristles at the side of the abutment rather than at the skin. Your focus should be to remove crust and debris around the base of the

abutment rather than brushing your skin. Your skin is delicate and should not be treated too harshly.

Remember, it is important to keep your soft brush clean too!

After cleaning, gently dry the area around the abutment with a clean towel or tissue. Do not allow your hair to wrap itself around the abutment. If you use a hairdryer do not keep the hot air on the abutment for too long.

Weekly care

If you use an abutment cover it will help keep the inside of the abutment clean however, the inside of the abutment should still be cleaned on a weekly basis. The soft brush can also be used for cleaning inside the abutment. Again, this is best done in conjunction with a bath or shower when soap and warm water can soften dirt making it easy to remove.

Skin debris or any other dirt in the abutment can prevent a smooth fitting of the abutment to your abutment cover and your Baha sound processor. This could potentially lead to damage of your Baha sound processor.

You will have routine check-ups initially every 4-6 months and later every 6-12 months.

However, should you experience any persistent irritation, soreness or any other signs of inflammation around the abutment, always contact your doctor.

Maintenance of your Baha Intenso™

The sound processor can be worn during waking hours apart from when you take a shower, or swim. Keep in mind that the sound processor is an electrical device that should be used with due care and attention. When engaging in physical activities, use common sense. If the activity is of such a nature that the sound processor risks getting damaged (for example contact sports)—remove it! If you choose to keep your sound processor on during physical activity we recommend securing it with the safety line (see below).

When not using your Baha Intenso™

When not in use ensure that the sound processor is switched off and stored carefully where it will not gather dust and dirt. If you will not be using the sound processor for any length of time, remove the battery. To maintain a good snap connection to the abutment, we recommend that you use an abutment cover when not wearing your sound processor.

Attachment of the safety line

If you keep the sound processor on when engaging in physical activities we advise you to attach a safety line to your sound processor so that it is not lost if it disconnects. The clip can then be fastened to clothing, for example, your shirt or jacket collar.

Cleaning your Baha


Ensure that your snap coupling is always visibly clean. Use a brush for cleaning purposes. (A regular toothbrush with hard bristles is suitable). Keep the plastic housing clean by regularly wiping with a soft, dry cloth.

Attention: For children who have been fitted with a Baha sound processor it is extremely important that the parent or caregiver of the child provides all necessary assistance and support in cleaning and maintenance, as detailed above.

Batteries

Zinc-air batteries type 675 should be used. Please note that these batteries are not rechargeable.

Battery replacement ► Open the battery compartment by sliding the cover downwards. Turn the sound processor over and, if the battery is stuck, shake the device gently. Place the new battery in with the + mark facing upwards. Press the battery down with your thumb and close the cover. If you are unable to slide the cover back easily, the battery may be incorrectly inserted. Please check that the + mark is facing upwards.

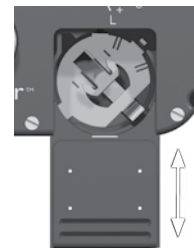Warning! The Baha sound processor contains small parts that may be a risk for small children. Batteries can be harmful if swallowed. Be sure to keep your batteries out of reach of small children or people of diminished mental capacity. In the event a battery is accidentally swallowed, seek medical attention at the nearest emergency centre immediately.

Battery life

Battery life begins decreasing as soon as the plastic strip is removed exposing it to air, even if the battery is not in use.

To maximise battery life, switch off the sound processor when not in use. Battery life varies depending on volume setting, the sound environment and the program setting. Battery life will decrease significantly when using the Baha sound processor with a connected FM receiver if the FM receiver uses the sound processor as its power source. Using powerful batteries is particularly important when high current consumption is common. The batteries included with the sound processor reflect our latest battery recommendation. Battery function can be related to a number of problems including: no sound, intermittent sound and crackling/buzzing. Thus, when encountering problems one of the first things one can do is to try changing the battery.

Note: Always follow local regulations for battery disposal.

General advice

- A sound processor will not restore normal hearing and will not prevent or improve a hearing impairment resulting from organic conditions.
- Infrequent use of a sound processor does not permit a user to attain full benefit from it.
- The use of a sound processor is only part of a hearing rehabilitation and may need to be supplemented by auditory training and instructions in lip-reading.
- For advice on "do's and don'ts" for your Baha sound processor see the section on Warnings and Precautions (p. 18).

Accessories

Several accessories are available to enhance your sound experience in different contexts. To ensure the smooth operation of your Baha, only use Cochlear original accessories or accessories approved by Cochlear.

Audio adapter

Allows direct input from computers, personal stereos, TVs, MP3s and other external equipment.

MicroLink Baha FM-receiver

Enables the use of most Phonak FM transmitters. For more information and orders contact your local Phonak distributor www.phonak.com

Telecoil unit

Improves sound quality and speech understanding when using telephones in homes where there are personal loop systems installed and in other buildings with loop facilities.

Direct audio input and microphone settings

The sound processor has a three-pin DAI (Direct Audio Input) port designed for connecting external equipment. As can be seen from the table below Program E allows sound from whatever device you have connected to the DAI only. In all the other programs, the microphone is still active and you will be able to hear sounds from your surrounding environment.

Attention: The marking on the plug must match up with the corresponding marking on the sound processor.

Safety instructions when connecting external equipment

Only connection cables and accessories that are supplied or approved by Cochlear should be used. Non-standard connection cables can cause damage to the sound processor and injury to you. If in doubt, consult your audiologist. Cochlear is not liable for any damages arising as a result from this or any other form of incorrect use of the product or its parts.

Warning! It is extremely dangerous and potentially fatal to connect any equipment that is plugged into a 110 and/or 220 V power supply without using the audio adapter.

Program	Microphone setting	DAI
1	on	on
2	on	on
E	off	on

Abutment cover

An abutment cover is included in the sound processor kit. The cover helps conceal and protect the abutment when the sound processor is not in place. This improves aesthetics and prevents dirt from collecting in the abutment, which helps maintain a good snap connection between the abutment and sound processor. Simply press the cover gently onto the abutment.

Share the Experience

The test rod provided can be used by others so that they can listen to your Baha sound processor themselves. Snap the sound processor onto the test rod and hold the rod against the skull bone behind the ear while keeping the ears plugged (take care not to allow the sound processor to touch anything except the test rod as this may cause feedback),

Warnings and precautions

- The sound processor is a digital, electrical, medical device designed for specific use. As such, due care and attention must be exercised by the user at all times.
- The sound processor is not waterproof! Never wear it in the bath or shower!
- Do not expose the sound processor to extreme temperatures etc. For example do not wear it in the sauna. The sound processor is designed to operate within the temperature ranges +5°C (+41°F) to +40°C (+104°F). In particular, battery performance deteriorates in temperatures under +5°C. The device should not be subject, at any time, to temperatures below -20°C (-4°F) or above +50°C (+122°F).
- This product is not suitable for use in the presence of flammable materials.
- If you are to undergo an MRI (Magnetic Resonance Imaging) procedure, the sound processor should be removed beforehand. The implant and abutment remain in place.
- Dispose of electronic items in accordance with your local regulations.
- Do not dispose of in fire.

Troubleshooter

If you have encountered any problems with your Baha sound processor:

- Check that plastic housing is not touching anything (such as a hat, eyeglasses or similar items).
- Check the volume control.

- Check the program setting
- Try changing the battery.

If the problem persists, contact your audiologist to get it fixed. Do not try fixing the sound processor yourself.

Service, repairs & warranty

Service and repairs

If there is any fault with your Baha sound processor, contact your ear clinic or hearing aid centre and they will help you to get it serviced.

Telephone number of hearing centre: _____

Contact person: _____

Telephone number of ear clinic: _____

Contact person: _____

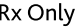
Warranty

The warranty does not cover defects or damage arising from, associated with, or related to the use of this Product with any non-Cochlear processing unit and/or any non-Cochlear implant. See "Cochlear Baha Global Limited Warranty card" for more details.

DISCLAIMER: Cochlear reserves the right to make any necessary changes and improvements to the product.

Key to symbols

Throughout the manual you will see small symbols. Please refer to the list below for explanations:

Symbol	Explanation	Symbol	Explanation
	"Consult instructions for use" or "Consult operating instructions"		Manufacturer
	"Caution, consult accompanying documents" or "Attention, see instructions for use"		Date of manufacture
	Audible signal		Batch Code
	CE-mark		Serial Number
	Temperature limitation		Caution. Federal law (USA) restricts this device to sale by or on the order of a licensed medical practioner
	Interference risk		

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Adresse :

Customer Service - Cochlear Americas

13059 E. Peakview Avenue, Centennial, CO 80111, USA

Toll free (North America) 1800 523 5798

Tel: +1 303 790 9010 Fax: +1 303 792 9025

Email: customer@cochlear.com

Customer Service - Cochlear Europe

6 Dashwood Lang Road, Bourne Business Park, Addlestone, Surrey KT15 2HJ, UK

Tel: +44 1932 87 1500 Fax: +44 1932 87 1526

Email: info@cochlear.co.uk

Customer Service - Cochlear Asia Pacific

1 University Avenue, Macquarie University, NSW 2109, Australia

Toll free (Australia) 1800 620 929

Toll free (New Zealand) 0800 444 819

Tel: +61 2 9428 6555 Fax: +61 2 9428 6352 or Toll free 1800 005 215

Email: customerservice@cochlear.com.au

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Instructions for use in the USA

Caution: Federal law (USA) restricts this device to sale by or on the order of a medical practitioner.

Warning to hearing aid dispensers

A prospective user of Baha® must always have a medical evaluation by a licensed physician (preferably an ear specialist) and have an implant placed in the skull bone before a Baha sound processor can be dispensed.

Important notice for prospective users of a Baha® sound processor

- A prospective user of Baha must always have a thorough medical evaluation by a licensed physician (preferably a physician who specialises in diseases of the ear, otolaryngologist, otologist, or otorhinolaryngologist) and have an implant placed in the skull bone before a Baha sound processor is purchased.
- The purpose of the medical evaluation is to ensure that all medically treatable conditions that may affect hearing are identified and treated before the sound processor is purchased. Following the medical evaluation, the physician will give you a written statement that states that your hearing loss has been medically evaluated and that you may be considered a candidate for Baha. The physician will refer you to an audiologist or a hearing aid dispenser, as appropriate, for an audiological evaluation.
- Federal law restricts the sale of Baha sound processors to those individuals who have obtained a medical evaluation from a licensed physician. Federal law permits a fully informed adult to sign a waiver statement declining the medical evaluation for religious or personal beliefs that preclude consultation with a physician. The exercise of such a waiver is not in your best health interest and its use is strongly discouraged.

A word of caution about batteries

Batteries can be harmful if swallowed. Be sure to always keep your batteries out of reach of small children or people of diminished capacity. In the event a battery is accidentally swallowed, seek prompt medical attention at the nearest emergency centre, or call the **National Button Battery Hotline collect at (202) 625 - 3333**

Manufacturer:

Cochlear Bone Anchored Solutions AB Konstruktionsvägen 14, SE-435 33 Mölnlycke, Sweden Tel: +46 31 792 44 00 Fax: +46 31 792 46 95

Regional offices:

Cochlear Americas 13059 E. Peakview Avenue, Centennial, CO 80111, USA Tel: +1 303 790 9010 Fax: +1 303 792 9025

Cochlear AG European Headquarters, Peter Merian-Weg 4, CH-4052 Basel, Switzerland Tel: +41 61 205 0404 Fax: +41 61 205 0405

Cochlear Ltd (ABN 96 002 618 073) 1 University Avenue, Macquarie University, NSW 2109, Australia Tel: +61 2 9428 6555 Fax: +61 2 9428 6352

Local offices:

Cochlear Benelux NV Schaliënhoedreef 20 i, B-2800 Mechelen, Belgium Tel: +32 15 79 55 77 Fax: +32 15 79 55 70

Cochlear Canada Inc 2500-120 Adelaide Street West, Toronto, ON M5H 1T1, Canada Tel: +1 416 972 5082 Fax: +1 416 972 5083

Cochlear (HK) Ltd Beijing Representative Office Unit 2208-2212, Tower B, No.9 Gemdale Building, No.91 Jianguo Road, Chaoyang District, Beijing, PR China 100022
Tel: +86 10 5909 7800 Fax: +86 10 5909 7900

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Cochlear Europe Ltd 6 Dashwood Lang Road, Bourne Business Park, Addlestone, Surrey KT15 2HJ, United Kingdom Tel: +44 1932 87 1500 Fax: +44 1932 87 1526

www.cochlear.com

Baha, Baha Divino, Baha Intenso and Vistafix are trademarks of Cochlear Bone Anchored Solutions. Cochlear and the elliptical logo are trademarks of Cochlear Limited. © Cochlear Bone Anchored Solutions 2010. All rights reserved.
MAR11. Printed in Sweden.

E80773E

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Cochlear Nordic AB Konstruktionsvägen 14, SE-435 33 Mölnlycke, Sweden Tel: +46 31 335 14 61 Fax: +46 31 335 14 60